


BiTrac Shield™

innovative technology for Non-Invasive Ventilation

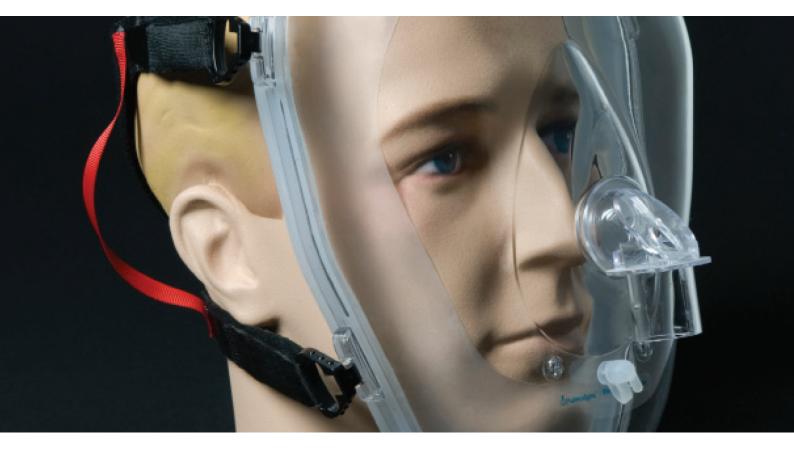

Pulmodyne[®] BiTrac Shield[™] is ideal for use in Non-Invasive Ventilation (NIV), CPAP and pressure support Bilevel and Ventilator Applications.

The BiTrac Shield provides a perfect solution for patients who can not tolerate full face or nasal masks and is ideal for patients requiring immediate ventilation.



Complete line of masks for immediate non-invasive ventilation

The **BiTrac Shield**[™] is an ideal solution for:

- Cross contamination concerns
- Cleaning and disinfection issues
- Claustrophobic Patients
- Patients who have developed nasal skin breakdown
- Patients with a facial deformities
- Edentulous Patients
- BiTrac Shield[™] for CPAP, Bilevel and Ventilator Applications
- Silicone cushion for immediate low pressure seal
- Prevents nasal skin breakdown
- Unobstructed view for claustrophobic patients
- One size fits most adults, simplifies mask fitting
- MRI compatible
- DISPOSABLE Single Patient Use
- Latex Free
- Made in USA

BiTrac Shield™

BiTrac Shield[™] allows Anti-Asphyxia Elbow for patient safety unobstructed view minimizing during Bilevel Applications (BiTrac SE claustrophobia. Shield has a Standard Elbow for ventilator application). Convenient, pre-attached Head Strap. BiTrac SE Shield[™] Manometer port for Head Strap with connecting proximal quick and easy pressure line or for attachment clips. monitoring CO₂, Two exhaust ports provide continuous leak (BiTrac SE Shield does not have exhaust ports).

Ordering information

BiTrac Shield [™] product information		Part numbers
BiTrac Shield [™] (Bilevel and CPAP Application) with Anti-Asphyxia Elbow 22mm Female and Shield Head Strap		
Adult	10 per case	313-9050
BiTrac SE Shield [®] (Ventilator Application) with Standard Elbow 22mm Female and Shield Head Strap		
Adult	10 per case	313-9051
Circuits product information		Part numbers
BiTrac NIV Circuit (Bilevel/CPAP Applications)		
72" Smooth ID Tubing with swivel connector, with exhalation port, cap plug; includes optional 24" tubing with "T"	10 per case	313-9074
BiTrac NIV Circuit (Bilevel/CPAP Applications)		
72" Smooth ID Tubing with attached filter, swivel connector with exhalaiton port, cap plug; includes optional 24" tubing with "T"	10 per case	313-9075
ValueTrac NIV Circuit (Bilevel/CPAP Applications)		
72" Corrugated Tubing and swivel connector with exhalation port	10 per case	313-9078