

# VacuAide 7325 Portable Suction Unit

For effective airways clearance



**Li-ion battery provides a fast charging time and a long battery life.**

**Very lightweight at only 1.7kg including the battery.**

**The VacuAide 7325 Portable Suction Unit has been designed to perform effective airways suctioning in a comfortable, safe and hygienic way.**

Featuring a lightweight Li-ion battery with a quick charging time and long battery life for improved performance. The VacuAide 7325 is suitable for both a clinical and homecare setting.




The 7325P is supplied with its own carry bag designed to allow discreet portable use of the device.



For convenience and hygiene an integrated tube holder allows the user to place the suction tube after each suction procedure.

**drive**

 **DeVilbiss**  
HEALTHCARE

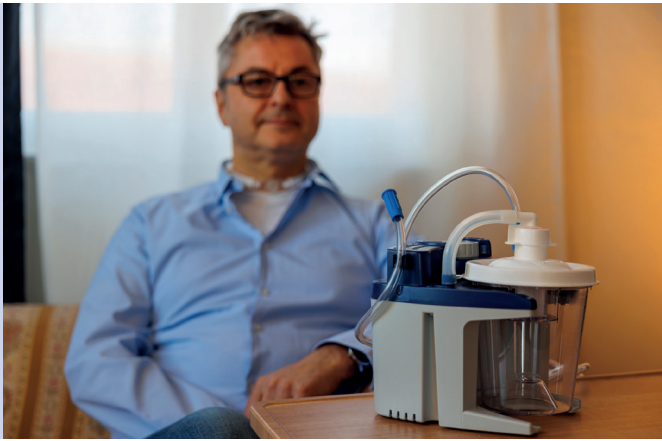

## Features

- Features an internal Li-ion battery (7325P) with a fast charge time of 3 hours, a runtime of 45 minutes and a 1 year warranty on the battery
- Very lightweight at 1.7kg including the battery for easy transportation
- A suction tube holder is integrated into the 7325 for convenience and hygiene
- The vacuum gauge is shock protected and can be adjusted from 50 - 550mmHg
- Easy to operate and maintenance free
- Can be operated discreetly inside the personal carry case
- The personal carry case includes a viewing panel so that you can visually check the secretion level at any time. It is padded and is designed as a shoulder bag, with internal compartments for accessories
- Complete with 100-240 V power supply and 12 V DC car adaptor
- Collection canister is easy to install and remove
- The VacuAide 7325 comes complete with an 800ml collection canister
- A reusable autoclavable 1200ml canister and the Flovac canister and disposable liner system are available as accessories
- LED control lights are displayed on the top of the device for quickly checking the units charge status
- High suction power – up to 27 l/min free flow
- To prevent liquids entering the suction unit there are two safety features; a hydrophobic filter and a mechanical overflow protection built into the canister lid
- Shares many of the consumables with the existing VacuAide suction devices, including collection canisters, filters and patient tubing
- A mains powered version is also available (D model). Please note this doesn't include the internal battery, 12 V car adaptor or carry case

**If you would like to view this product, please contact Air Liquide Healthcare.**

For more information please contact us at:

Air Liquide Healthcare, Stratus House, Unit 1 IDA College Park, Blanchardstown Road North, Blanchardstown, Dublin 15, D15 PEC4.

Tel: +353 (0)1 809 1800 Freephone: 1800 24 02 02 (ROI only) 0800 328 5875 (NI only)

Fax: +353 (0)1 829 3966 Email: [healthie@airliquide.ie](mailto:healthie@airliquide.ie) Website: [www.airliquidehealthcare.ie](http://www.airliquidehealthcare.ie)

## Product Codes

Description	Part number
VacuAide 7325 Portable Suction Unit (with internal rechargeable battery)	<b>7325P-U</b>
VacuAide 7325 Mains Powered Suction Unit (without internal battery, 12V car adaptor or carry case)	<b>7325D-U</b>

Product specifications	
Height	17.5cm
Width	17.8cm
Depth	18cm
Operating Temperature Range:	0°C - 40°C
Product Weight 7325P-U	1.7kg
Product Weight 7325D-U	1.5kg
Air Flow	27 l/min
Vacuum Range	50 – 550mmHg
Charging time (for 7325P)	Up to 3 hours
Battery running time (for 7325P)	Up to 45 mins
Device Warranty	2 Years
Battery Warranty	1 Year
Electrical Requirements	100-240V AC, 50/60 Hz; 12V DC

